

## Optional display for the BlueInOne


Measurement data for CO<sub>2</sub> and O<sub>2</sub> always in sight

For all the BlueInOne gas analyzers there is a convenient additional display available now. It shows the current measurement data of the CO<sub>2</sub> and O<sub>2</sub> concentration. So you have the present measurement always directly in view. The display is simply connected via an up to five meters long data cable to the analogue 4-20 mA output of the BlueInOne analyzer. Magnets at the backside of the housing of the display make the mounting easy and convenient. The display also has got a 4-20 mA port, which passes the analog signal from BlueInOne. Thus, no ports are blocked by the display. In addition, the display has a practical button to perform a 1-point calibration. Thus, the calibration can also be started very simply, if the BlueInOne is not easy to reach. The display is supplied via the power from the 4-20 mA output and requires no external power supply and no batteries. There are versions for each BlueInOne; also devices that are already in use can be upgraded.



### Advantages:

- Convenient: the current measured values directly on site in view
- Easy installation: quick control without software
- Flexible to install, easy mounting
- Time savings: start the one-point calibration easily even if the BlueInOne is hard to reach
- 4-20 mA port included, no ports are blocked
- Variants for all BlueInOne available

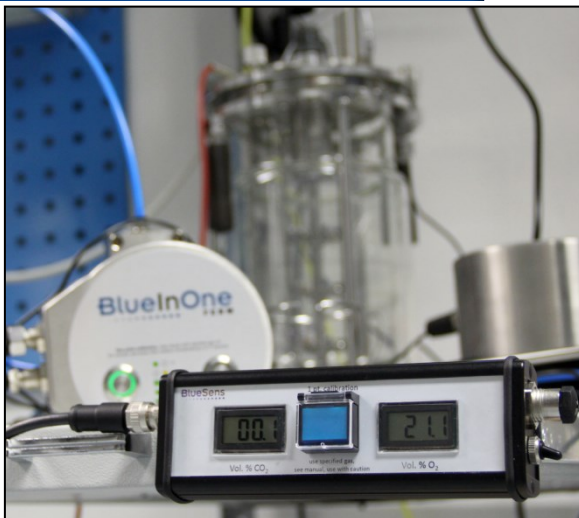

Keep the readings on site in view



Easy installation